# ROWENTA®

goodvibes<sup>®</sup>  
LIFE+

BR972x

Instruction manual - Notice d'utilisation - Gebruiksaanwijzing - Bedienungsanleitung - Manuale di istruzioni - Manual de instrucciones - Manual de instruções - Εγχειρίδιο οδηγιών - Kullanma kılavuzu - Instrukcja obsługi - Uživatelská příručka - Használati útmutató - Návod na použitie - Manual de instruții - Uputstvo za upotrebu - Priručnik s uputama - Ръководство синструкци - Руководство по эксплуатации - Priročnik znavodili - Посібник з експлуатації - Пайдалану нұсқаулығы Ръководство с инструкции - Priručnik s uputstvima - Manual de instrucciones - Kasutusjuhend - Használati útmutató - Manuale di istruzioni - Naudojimo instrukcija - Lietošanas rokasgrāmata - Manual de instruções - Navodila za uporabo - Uputstvo za upotrebu - Talimat kılavuzu

[www.rowenta.com](http://www.rowenta.com)

1820017014/01



EN  
FR  
NL  
DE  
RU  
UK  
RO  
BU  
PL  
HR  
CS  
SK  
KK  
EL  
BG  
BS  
ES  
ET  
HU  
IT  
LT  
LV  
PT  
SL  
SR  
TR

1

**EN** - Body composition scale / **FR** - Impédancemètre / **NL** - Impedantiemeter  
**DE** - Körperfettanalysewaage / **RU** - Весы с измерением жировой мышечной масс  
**UK** - Ваги - аналізатор складу тіла / **RO** - Cântar corporal  
**BU** - Персонална весна със скала за състава на тялото / **PL** - Waga z funkcją analizy składu ciała  
**HR** - Vaga za građu tijela / **CS** - Osobní diagnostická váha / **SK** - Váha s analýzou stavby tela  
**KK** - Дене құрамын өлшеу таразысы / **EL** - Συγάρια σύστασης σώματος  
**BG** - Персонална весна със скала за състава на тялото / **BS** - Vaga za kompoziciju tijela  
**ES** - Báscula de medición corporal / **ET** - Kehakoostise analüüsi kaal / **HU** - Testösszetétel-mérleg  
**IT** - Bilancia impedenziometrica / **LT** - Svarstyklės su kūno analizės funkcija  
**LV** - Ķermeņa uzbūves svāri / **PT** - Balança de composição corporal  
**SL** - Tehtnica za merjenje telesne sestave / **SR** - Vaga za merenje telesnog sastava  
**TR** - Vücut bileşimi tartisi

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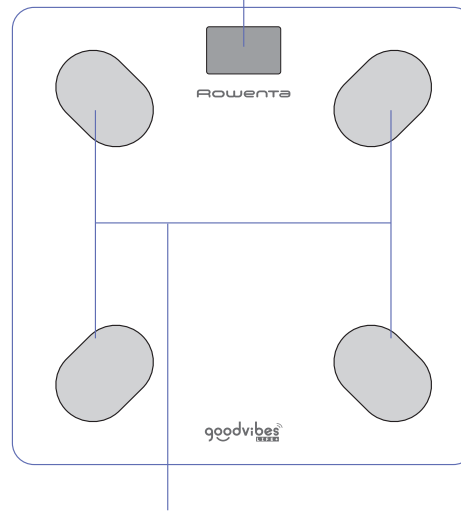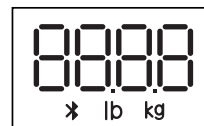**EN** Thanks for choosing TEFAL Goodvibes body composition Scale. This product employs BIA (bioelectrical impedance analysis) technology for your body composition data and Bluetooth technology to connect the scale to your smartphone. To have accurate and consistent measurement results, please use it at the same time (preferably twice a week, 15 min after waking up) and under same conditions. Please avoid measurements:

- During or immediately after strenuous exercises
- Immediately after a shower, bath or sauna
- With heavy alcohol consumption
- During pregnancy
- During fever

During the measurement, please stand straight and barefoot. Wearing socks or dirty feet can affect the measurement accuracy. The scale should not be used to diagnose or treat any medical condition. All data is for reference and can not be a substitute for advice from a licensed healthcare professional.

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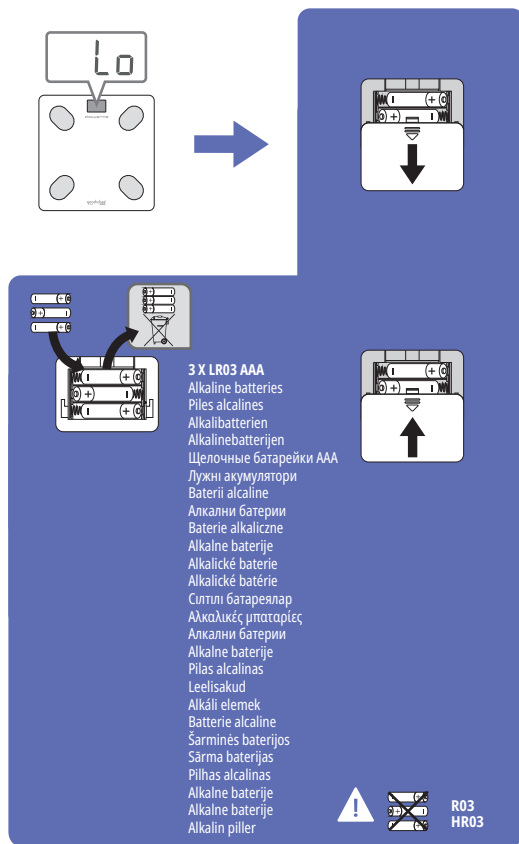
**EN WARNING**

- Do not attempt to recharge non-rechargeable batteries.
- Do not mix new and used batteries or different types of batteries.
- Ensure that the polarity is correct when inserting the batteries.
- Empty batteries must be removed from the device and disposed of in a safe manner.
- Remove batteries from a device if it is not intended to be used for a long period.
- Do not short-circuit the power supply terminal.

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