Tefal®

EN

FR

NL

DE

RU

goodvibes[®]
HEALTH

BM966x


Instruction manual – Notice d'utilisation – Gebruiksaanwijzing – Bedienungsanleitung – Manuale di istruzioni – Manual de instrucciones – Manual de instruções – Εγχειρίδιο οδηγιών – Kullanma kılavuzu – Instrukcja obsługi – Uživatelská příručka – Használati útmutató – Návod na použitie – Manual de instrucțiuni – Uputstvo za upotrebu – Uputeza uporabu – Ръководство с инструкции – Руководство по эксплуатации – Priročnik znavodili – Посібник з експлуатації

1

EN - Body composition scale/ **FR** - Impédancemètre/ **NL** - Impedantiemeter/
DE - Körperfettanalysewaage/ **RU** - Весы с измерением жировой мышечной масс

EN - Weighing unit modification can only be done in the app
FR - Le changement d'unités de poids n'est possible qu'à partir de l'application
NL - De weegeenheid kan alleen via de app worden gewijzigd
DE - Die Änderung der Wiegeeinheit kann nur über die App vorgenommen werden
RU - Единицы измерения массы можно настроить только в приложении

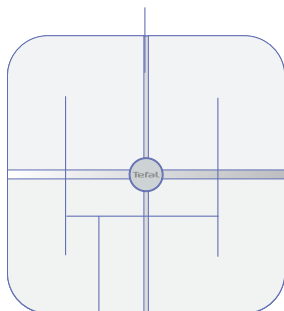
EN Electrodes

FR Electrodes

NL Electroden

DE Elektroden

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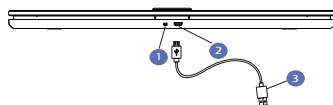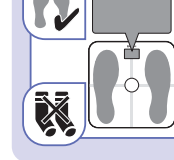
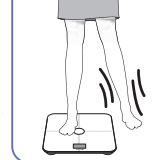
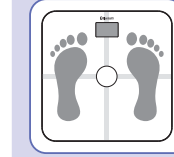
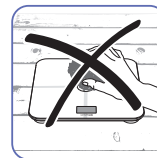
EN Thanks for choosing TEFAL Goodvibes body composition Scale. This product employs BIA (bioelectrical impedance analysis) technology for your body composition data and Bluetooth technology to connect the scale to your smartphone. To have accurate and consistent measurement results, please use it at the same time (preferably twice a week, 15 min after waking up) and under same conditions. Please avoid measurements:

- During or immediately after strenuous exercises
- Immediately after a shower, bath or sauna
- With heavy alcohol consumption
- During pregnancy
- During fever

During the measurement, please stand straight and barefoot. Wearing socks or dirty feet can affect the measurement accuracy. The scale should not be used to diagnose or treat any medical condition. All data is for reference and can not be a substitute for advice from a licensed healthcare professional.

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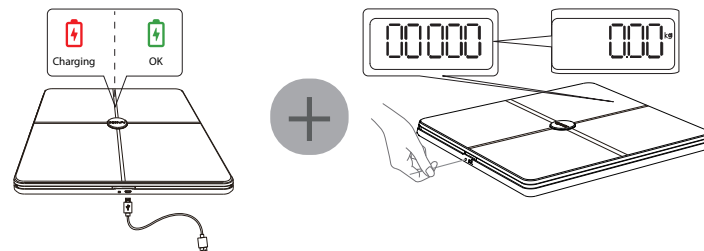
EN

Launch the App

Set the weighing units (kg / st.lb / st) and the height units (cm / inch)
Register your profile user

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
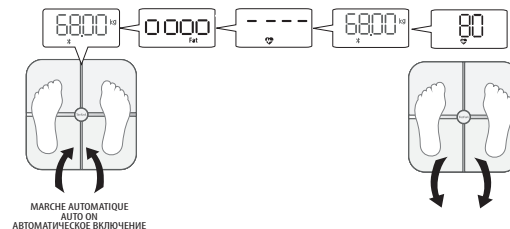
EN - Please wait until all the data is displayed in the App to step off.

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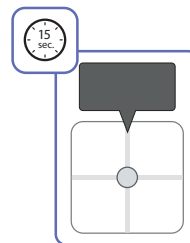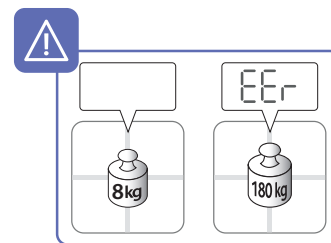
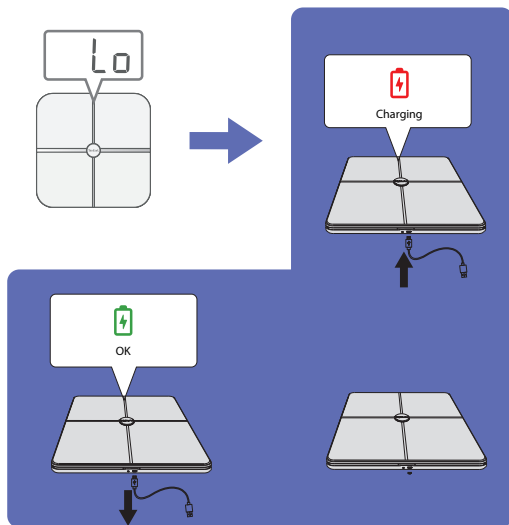
EN This appliance contains batteries that are only replaceable by skilled persons.

WARNING:

- Replacement of a battery with an incorrect type can defeat a safeguard.
- Disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery can result in an explosion.
- Leaving a battery in an extremely high temperature surrounding environment can result in an explosion or the leakage of flammable liquid or gas.
- A battery subjected to extremely low air pressure can result in an explosion or the leakage of flammable liquid or gas.
- For the purposes of recharging the battery, only use a detachable supply unit USB 5V 0,5A.

ENVIRONMENT PROTECTION FIRST !

- ① Your appliance contains valuable materials which can be recovered or recycled.
➡ Leave it at a local civic waste collection point.

RECYCLING THE BATTERIES: Before removing the Li-ION batteries, discharge them completely. Open the battery compartment using a tool, cut all the connections and then recycle the batteries in a safe manner.

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